

Implemented Samples Page 18/09/2026

- Added rowSelectionModel and dialogs after deletion. These are different features compared with Batches Page. Implemented these features also for Batches Pages.

<

Dashboard

Experiment

Acquire

Equipment

Data

Batches

Samples

Files

Analysis

Graphs

Export

Settings

Database

Users

Chat

User

Samples

Search samples

Advanced Search

☒ Id Range

from

1

to

☐ Batch Id

List of found samples

	ID ↓	Sample Name	Description	Batch ID
☒	13	Dev14	1000 and 2000 um in diameter diodes. ...	20
☐	10	Dev26	MBE30_26 on MgO Si11	10
☐	9	Dev06	MBE26_06 on c-Sap Si 10	10
☐	8	Dev20	300, 500, and 800 um in diameter diode	12
☐	7	Dev08	1000 and 2000 um in diameter diode	12

1 row selected

Rows per page: 5 1-5 of 11

Sample ID

13

Name*

Dev14

Description

1000 and 2000 um in diameter diodes. O5632 NW14.

Treatment

Standard treatment.

Properties

Batch ID

20

* = Required

NEW

ADD

UPDATE

DELETE

© 2026 ExpeGraph by Giang T. Dang.

Implemented Search batch button with advanced search, 18/08/2026 database in desktop work computer

<

Dashboard

Experiment

Acquire

Equipment

Data

Batches

Samples

Files

Analysis

Graphs

Export

Settings

Database

Users

Chat

Batches

Search batches

Advanced search

ID range

from

to

Fabrication date range

from

01/01/2024

to

01/12/2024

Project Id

Lab Id

List of found batches

	ID	Batch Name	Fabrication Date	Treatment
	1	IrOxNewSM	2024-12-01	Standard tr...
	10	MBEmesfets	2024-12-01	Standard tr...
	12	IrOxNewSMAftCoat	2024-12-01	Standard tr...

Rows per page:

5

1-3 of 3

<

>

Batch ID

Batch Name *

IrOxNewSM

Description

IrOx SBDs using new milled shadow masks

Fabrication Date *

01/12/2024

Treatment

Standard treatment: Ohmic 450 C, 3min, N2, furnace.
IrOx/Ir 40/40 nm

Project ID

Lab ID

* = Required

NEW

ADD

UPDATE

DELETE

© 2026 ExpeGraph by Giang T. Dang.

Implemented Search batch button with advanced search, 18/08/2026 database in desktop work computer

<

Dashboard

Experiment

Acquire

Equipment

Data

Batches

Samples

Files

Analysis

Graphs

Export

Settings

Database

Users

Chat

User

Batches

Search batches

Q

Advanced search

☒ ID range

from

1

to

20

☒ Fabrication date range

from

01/01/2024

to

18/12/2024

☐ Project Id

☐ Lab Id

List of found batches

☐

ID

Batch Name

Fabrication Date

Treatment

☐

1

IrOxNewSM

2024-12-01

Standard tr...

☐

10

MBEmesfets

2024-12-01

Standard tr...

☐

11

NegLOBefCoat

2024-12-05

Standard tr...

☐

12

IrOxNewSMAftCoat

2024-12-01

Standard tr...

☐

13

NegLOAftCoat

2024-12-05

Standard tr...

☐

14

NewSMnewGO

2024-12-06

Standard tr...

☐

15

Thicklr

2024-12-16

Standard tr...

Rows per page: 50

1-7 of 7

Batch ID

Batch Name *

Description

Fabrication Date *

Treatment

Project ID

Lab ID

* = Required

NEW

ADD

UPDATE

DELETE

Added new
Grid views for
results of
searched
batches
25/07/2026
Database in
Laptop
computer

<

Dashboard

Experiment

Acquire

Equipment

Data

Batches

Samples

Files

Analysis

Graphs

Export

Settings

Database

Users

Chat

Batches

Search batches

Q

> Advanced search

List of batches

☐	ID	Batch Name	Fabrication Date	Treatment	
☐	1	batch 1	01/12/2024	treatment A	
☐	2	batch 2	02/01/2025	treatment A	
☐	3	batch 3	02/02/2025	treatment A	
☐	4	batch 4	02/03/2025	treatment A	

Rows per page: 5 1-5 of 8 < >

Batch ID

Batch Name *

Description

Fabrication Date

Treatment

Project ID

Lab ID

* = Required

NEW

ADD

UPDATE

DELETE

User

© 2026 ExpeGraph by Giang T. Dang.

In-Confidence

Change
Icons
10/07/2026

Dashboard

Experiment

Acquire

Equipment

Data

Batches

Samples

Files

Analysis

Graphs

Export

Settings

Database

Users

Chat

Acquire

SETUP

RUN

Connected Equipment

Select equipment

Measurement Range

V_{start}

0

V_{stop}

1

V_{step}

0.1

+

Display Setup

X Axis

Mode

Auto

Scale

Linear

Y Axis

Mode

Auto

Scale

Linear

Stop
Button
5/5/2026

Dashboard

Experiment

Acquire

Equipment

Data

Batches

Samples

Files

Analysis

Graphs

Export

Settings

Database

Users

Chat

User ▾

Acquire


SETUP

RUN

RUN

STOP

I-V Curve

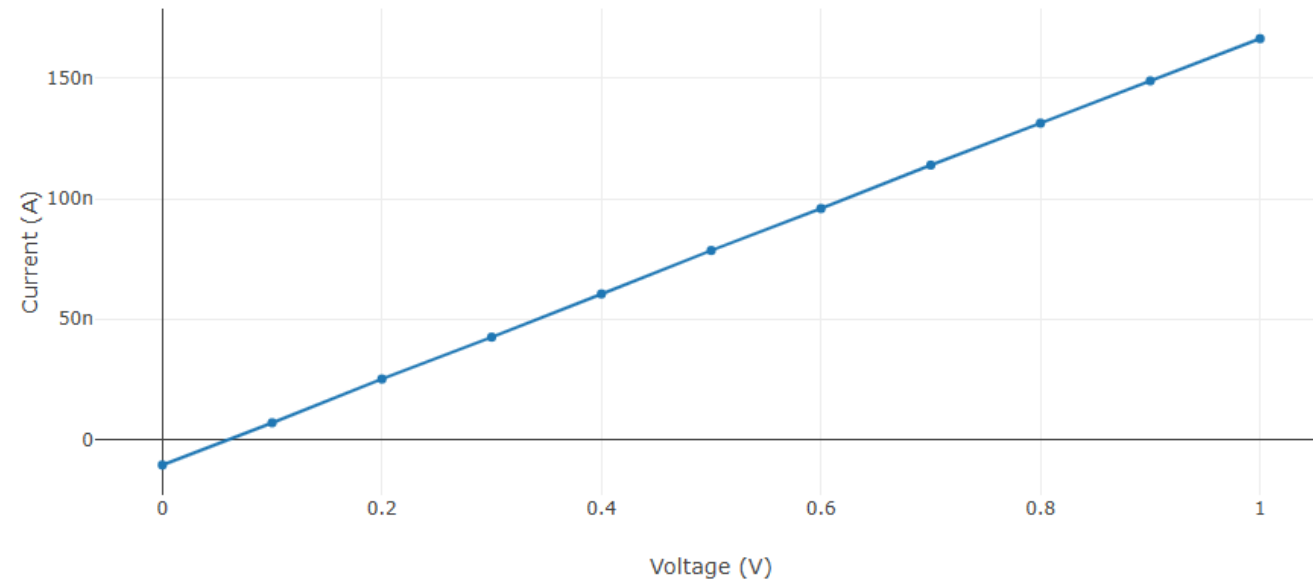

SAVE DATA TO .CSV

Save Data
To .CSV
Button
4/5/2026

Dashboard

Experiment

Acquire

Equipment

Data

Batches

Samples

Files

Analysis

Graphs

Export

Settings

Database

Users

Chat

User

Acquire


SETUP

RUN

RUN



I-V Curve

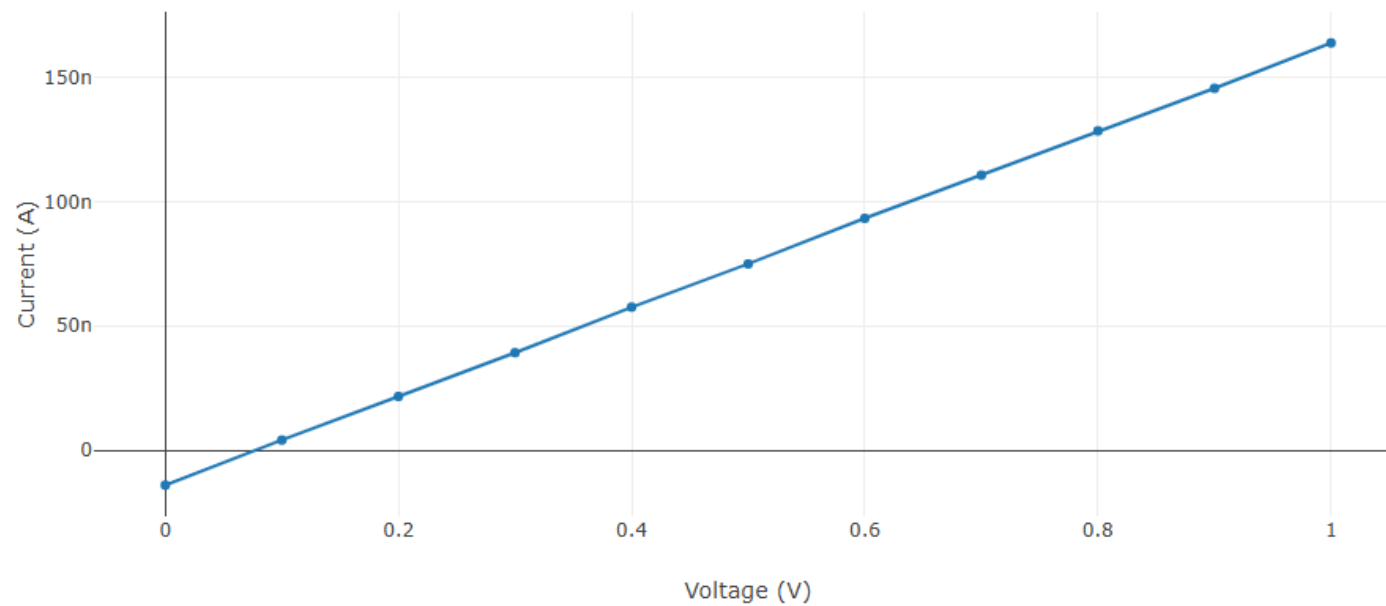

SAVE DATA TO .CSV

Dashboard

Experiment

► Acquire

⊗ Equipment

Data

■ Batches

◆ Samples

□ Files

Analysis

~ Graphs

↗ Export

Settings

○ Database

Ⓒ Users

 Chat

Acquire

SETUP

RUN

Connected Equipment

Keithley 6487 PICOAMMETER / VOLTAGE SOURCE ▾

Measurement Range

 V_{start}

0

 V_{stop}

1

 V_{step}

0.1

+

 V_{start}

0

 V_{stop}

-1

 V_{step}

0.1

-

 V_{start}

0

 V_{stop}

2

 V_{step}

0.1

-

 V_{start}

0

 V_{stop}

-2

 V_{step}

0.1

-

Display Setup

X Axis

Mode

Auto ▾

Scale

Linear ▾

Y Axis

Mode

Auto ▾

Scale

Linear ▾

Run
Button
18/4/2026

Dashboard

Experiment

Acquire

Equipment

Data

Batches

Samples

Files

Analysis

Graphs

Export

Settings

Database

Users


Chat

Acquire

SETUP

RUN

RUN

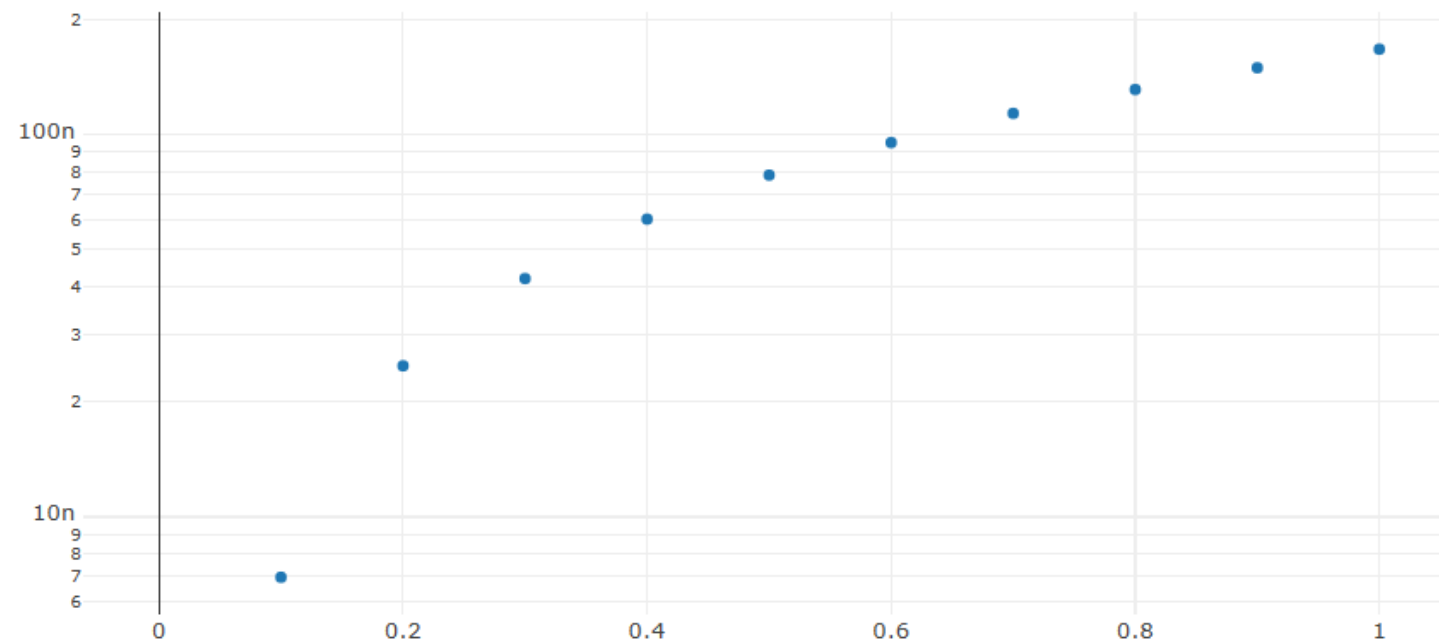

 Dashboard

Experiment

► Acquire

 Equipment

Data

 Batches

◆ Samples

□ Files

Analysis

~ Graphs

↗ Export

Settings

○ Database

◎ Users

 Chat

Equipment

All

New

Keithley 2461 SourceMeter®

Keithley 2657A HIGH POWER SYSTEM SourceMeter®

Keithley 2400 SourceMeter

Keithley 6487 PICOAMMETER / VOLTAGE SOURCE

Agilent B1500A Semiconductor Device Parameter Analyzer

CONNECT EQUIPMENT

DISCONNECT EQUIPMENT

Id

4

Name

Keithley 6487 PICOAMMETER / VO

Manufacturer

Tektronix

Model

6487

Serial Number

Purchase Year

Calibration Due

dd/mm/yyyy



Location

Connecting String

GPIB0::22::INSTR

Notes

ADD

UPDATE

DELETE